


ELAX Heating/Cooling  
Eye Massager

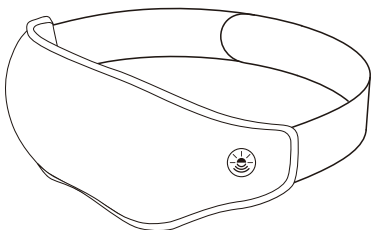

**USER MANUAL**

Please keep and read  
the manual before using.

# ABOUT

The ultra-soft and light-weighted eye mask fits perfectly to your eyes. It provides warm compression as well as high vibrating massage, helping to accelerate the blood circulation. Perfect for eye muscle fatigue.


Using with gel ice pack, it will help to reduce appearance of puffiness and relax the eyes.

Foldable & portable, wireless using, enjoy eye massage anytime and anywhere.

## SPECIFICATIONS

Product Name	Heating/Cooling Eye Massager
Model	ELAX
Operating voltage	3.7V ==
Working power	6W
Charging voltage	5V==
Charging current	2A
Heat temperature	45 ± 5°C
Batery capacity	600mAh
Product Size	240x95x75mm
Host weight	77.5g

## PART NAME

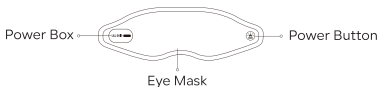

# HOW TO USE

1. Put the power box magnetically attached to the eye mask, after correct installation, the power indicator light will be on.



2. Fit the eye mask on eyes, and long press the power button on the right side of the eye mask for 1 second, the eye mask will be turned on and there will be a long vibration to indicate.



3. Power on the default massage mode (heat + vibration), short press the power button to switch between massage mode (heat+vibration) and heat mode (heat only), there will be 2 short vibrations to indicate.




4. If you need cold compress, please put the gel ice pack into the fridge for 10 minutes, then place it inside the eye mask to use without turning on the machine.



Caution: Automatic shutdown will happen after 18 minutes in each mode. If you need to stop the massage, you can press and hold the power button for 1s to shut down the mask. There will be a long vibration to indicate.

# HOW TO CHARGE

1. The power box is divided into four compartments. The lights will be on when the power box is magnetically attached to the eye mask, different light state represents to different power range. When the battery power is too low, the light will keep flashing.

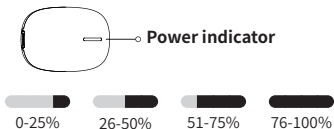

2. Before charging, please remove the power box from eye mask and charge it by the USB charging wire. Charging time is about 45 min. When charging, the blue indicator light will be flashing, and the light will stop flashing after fully charged. After charging finished, please unplug the charging cable in time.

## Failure Analysis

Failure	Possible Cause	Solution
Power Indicator light is not on	Poor contact between battery and mask Low Battery	Connect the battery well Charge the battery
Not working	Low battery Automatically off after 18 min of using	Charge the battery Re-start it
Cannot be charged	Poor contact between USB charging wire and power box Your adapter not working well	Connect the charging wire well Change another adapter

# WARNINGS

**To avoid any damage to person, please follow user instructions details outlined below.**

Consult your doctor before use if you are under medical treatment or if you have any of the following conditions:

- 1.Pregnant women.
- 2.People who have a body temperature of 38 °C or higher.
- 3.People who have lost sense of their skin temperature.
- 4.People who have suffered from low temperature burns.
- 5.People with abnormal body sensation.

People with the following conditions should not use this product.

- 1.People who are prohibited from exercise by medical advice.
- 2.Such as: thrombosis, severe aneurysms, acute venous aneurysms, all kinds of dermatitis, all kinds of skin infections. It is prohibited because it can easily cause deterioration.
- 3.People who cannot express their will or cannot operate the product by themselves.
- 4.Children under 14 years old cannot use this product.
- 5.If you feel dizziness, vomiting, throbbing or other physical reactions, please stop using it immediately.

# CAUTIONS

- 1.Do not use near the fireplace, or electric heaters and other warming appliances.
- 2.Do not use in a humid environment.
- 3.Do not use the product in the bathtub, shower, washbasin, swimming pool.
- 4.Do not place heavy objects on the product after use to prevent damage to the product.
- 5.Do not disassemble the product.
- 6.Avoid strong sunlight and corrosive gas and liquid erosion.
- 7.Keep this product out of reach of children.
- 8.The eye mask is NOT washable. Please use a clean wet cloth to wipe the surface of the mask.
- 9.Do not allow water to enter the product interior.
- 10.If it is not used for a long time, clean it up and put it in a storage bag.

Japan L&L Skin Tech Co., Ltd.  
Address: 3F Kabutocho Daiichi Heiwa Building 5-1  
Nihonbashi Kabuto-cho Chuo-ku, Tokyo, Japan 103-0026  
[www.llskinshop.eu](http://www.llskinshop.eu)



RoHS

